

OXYGEN HOMEFILL GUIDE

STEP 1

- Connect HomeFill (top unit) to oxygen concentrator (bottom unit) using curly supply tubing, and make sure that the HomeFill power cord and concentrator power cord are plugged in to the wall.



STEP 2

- Turn on oxygen concentrator (bottom unit)




STEP 3

- Set liter flow to patient's prescribed setting
- Look for green light to proceed
- *Setting cannot be above 2.5 for 5 L units or above 6 for 10 L units



STEP 4

- Remove cap from oxygen port

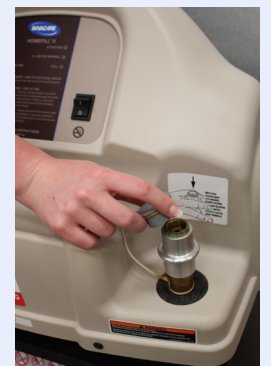


STEP 5

- Check tank to ensure cylinder valve is closed by turning clockwise.



STEP 6

- Before coupling the tank with the HomeFill, push down on the coupler collar until the GREEN dots are visible.

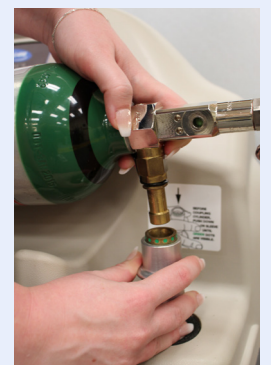

STEP 7

- Insert oxygen tank firmly, making sure it “clicks” into place.



STEP 8

- Turn on the HomeFill (top unit).



STEP 9

- The HomeFill will indicate a “BELOW NORMAL” or yellow light for 3 minutes prior to the HomeFill starting to fill the cylinder.



STEP 10

- When HomeFill is ready, it will switch to the “FILLING” or green light while oxygen cylinder is being filled. Filling could take several hours depending on the cylinder size.



STEP 11

- When oxygen cylinder is completely filled, the green “FILLING” light will go out and the “FULL” light will come on.



STEP 12

- Turn OFF the HomeFill.



STEP 13

- Remove the oxygen cylinder by pushing down on the HomeFill silver coupler collar and at the same time pull upward on the oxygen cylinder to unlock the connection and release it.



STEP 14

- If the concentrator is not being used, turn it off.



If a cylinder does not seem to fill properly, please start at the beginning of this guide and carefully repeat each step.

IMPORTANT NOTE – If the HomeFill “Attention” green light is on, contact Equipped for Life for assistance at **(301) 714-0200**